## ENGLISH

### Introduction

Read this user manual carefully before you use the straps for the BlueTouch Pain Relief Patch and save it for future reference.

### Intended use

The BlueTouch Pain Relief Patch is intended to be used on the upper or lower back and can be positioned and adjusted on the treatment area by means of a fabric strap.

### Important

- Only use the BlueTouch Pain Relief Patch in one of the Philips straps that has been specially designed for it. The straps are not supplied with the patch and have to be obtained separately.
- To prevent damage to the drum of the washing machine or to other laundry in the drum, wash the straps in a fine-mesh washing bag.

### Explanation of symbols

- This symbol means: Wash at a temperature of 30°max (Fig. 1).
- This symbol means: Do not bleach (Fig. 2).
- This symbol means: Do not tumble-dry (Fig. 3).
- This symbol means: Do not iron (Fig. 4).
- This symbol means: Do not dry clean (Fig. 5).
- This symbol means that the material has been tested for harmful substances according to Oeko-Tex® standard 100. (Fig. 6)
- This symbol means that this strap is to be used on the lower back. (Fig. 7)
- This symbol means that this strap is to be used on the upper back. (Fig. 8)

### Using your BlueTouch straps

The preferred method for fastening the BlueTouch Pain Relief Patch to the body is by means of one of the fabric straps. There is a special strap for treatment of the upper back and a special strap for treatment of the lower back.

### Positioning the lower back strap

To place the BlueTouch Pain Relief Patch in the strap, follow the steps below.

- 1 Put the lower back strap and the BlueTouch Pain Relief Patch in front of you on a table. Slide the left side of the BlueTouch Pain Relief Patch behind the elastic edge of the opening in the strap. (Fig. 9)
- 2 Then slide the other end of the BlueTouch Pain Relief Patch behind the opposite elastic edge of the opening so that the whole patch is inside the lower back strap. Check to make sure that the LEDs are not covered by the material of the strap. If necessary, shift the BlueTouch Pain Relief Patch a little (Fig. 10).
- 3 Remove clothing from the area to be treated.

*Note: Make sure the skin in the area to be treated is clean, entirely dry and free from creams, lotions and other cosmetics.*

- 4 Wrap the lower back strap with the BlueTouch Pain Relief Patch in it around your waist with the opening of the strap facing your back (Fig. 11).
- 5 Fasten it with the Velcro strips on both ends of the lower back strap. (Fig. 12)

*Note: Make sure the BlueTouch Pain Relief Patch in the strap is located on the area to be treated and the LED area is in contact with the skin.*

### Positioning the upper back strap

- 1 To put the BlueTouch Pain Relief Patch in the upper back strap, follow the same procedure as described in steps 1 and 2 in section 'Positioning the lower back strap' above.
- 2 Remove clothing from the area to be treated.
- 3 Put the upper back strap on your back with the opening of the strap facing your back (Fig. 13).

*Note: Make sure the BlueTouch Pain Relief Patch in the strap is located on the area to be treated and the LED area is in contact with the skin (Fig. 14).*

- 4 Fasten the waistband with the Velcro fasteners and fasten the shoulder straps to the two parts of waistband with the Velcro fasteners (Fig. 15).

*Note: You can also fasten the shoulder straps crosswise. This may provide a better fit in some situations (Fig. 16).*

### Cleaning

Always remove the BlueTouch Pain Relief Patch from the strap before you put the strap in the washing machine to clean it.

- The fabric straps can be washed in the washing machine (Fig. 17).

Caution: To prevent damage to the drum of the washing machine or to other laundry in the drum, wash the straps in a fine-mesh washing bag.

- Use a gentle cycle and a temperature of 30°C maximum. Dry the straps on the air and not in the tumble dryer. Do not iron the straps.

### Replacement

You can purchase new straps under typenumbers:

- PRP-LBS-6123: Lower back strap S/M
- PRP-LBS-6126: Lower back strap L/XL
- PRP-UBS-6113: Upper back strap S/M
- PRP-UBS-6116: Upper back strap L/XL

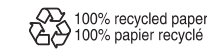
Register your product and get support at  
[www.philips.com/welcome](http://www.philips.com/welcome)

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